

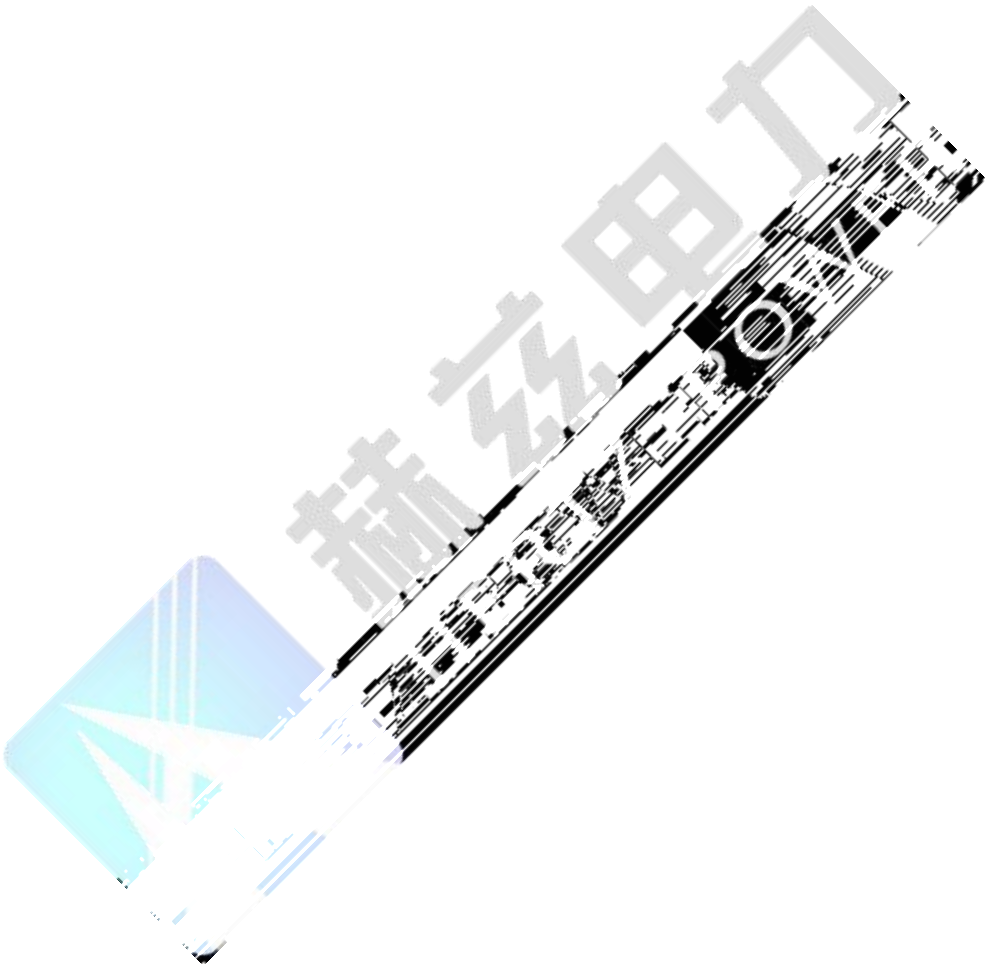


HZY-3S

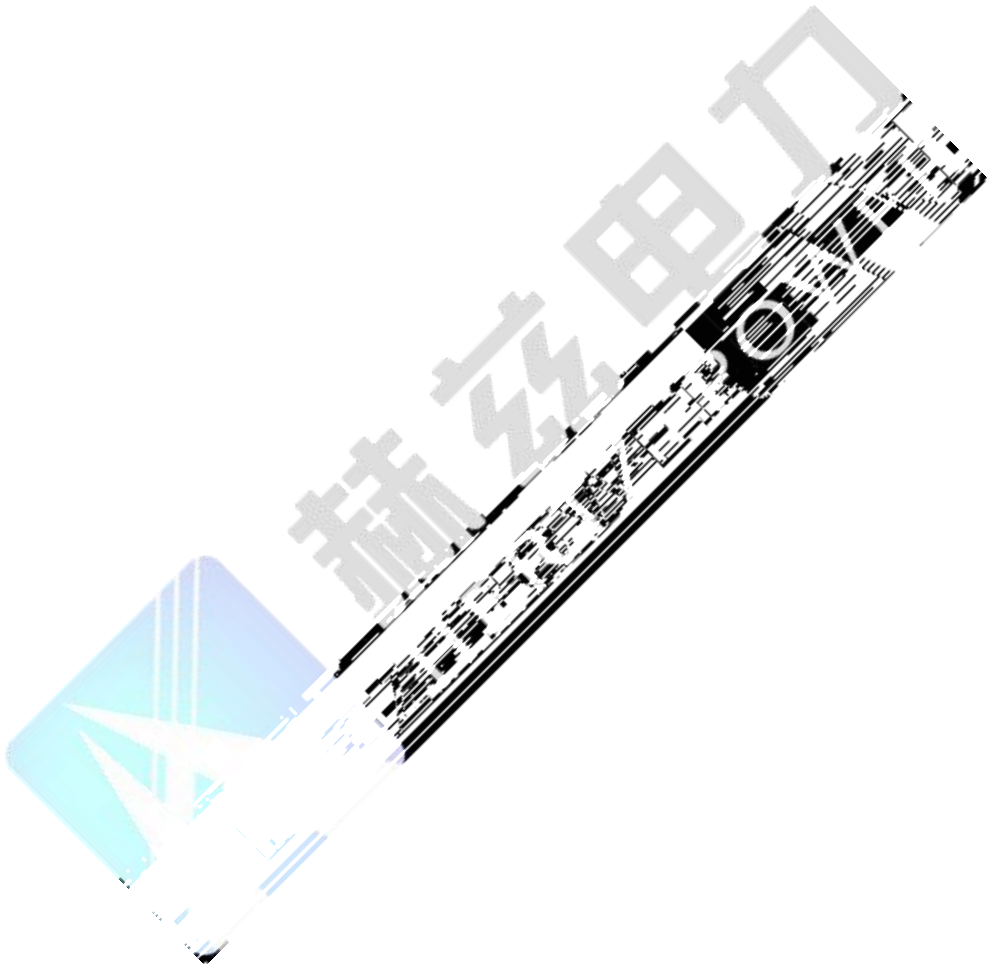

HZY-3S

1 3

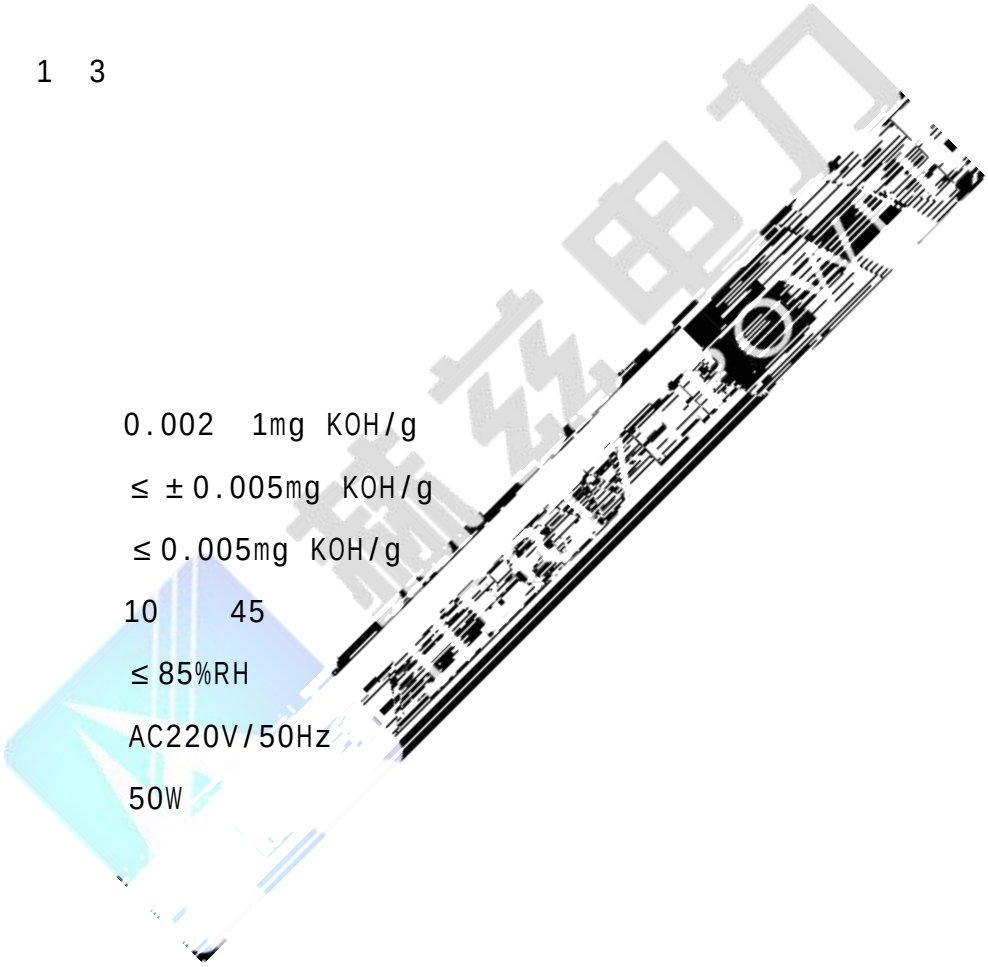

0.002 1mg KOH/g

≤ ±0.005mg KOH/g

≤0.005mg KOH/g

10 45

≤85%RH

AC220V/50Hz

50W



1

2 220V

220V

—KOH BTB

1 —KOH : KOH

—KOH 100ml

2 1000ml ( )

500ml

3 0.15mol/l

4 BTB BTB BTB


100ml 100ml

0.5ml 50ml BTB

5 KOH


50ml

6

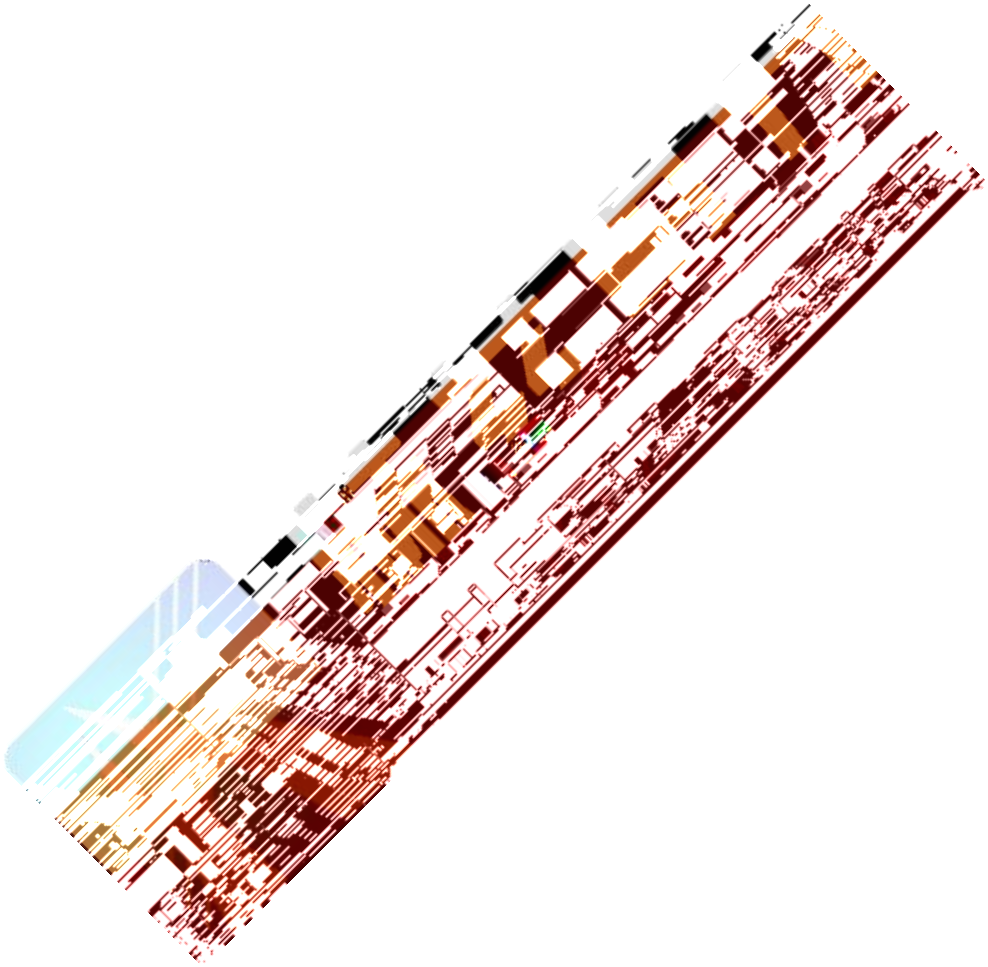





20

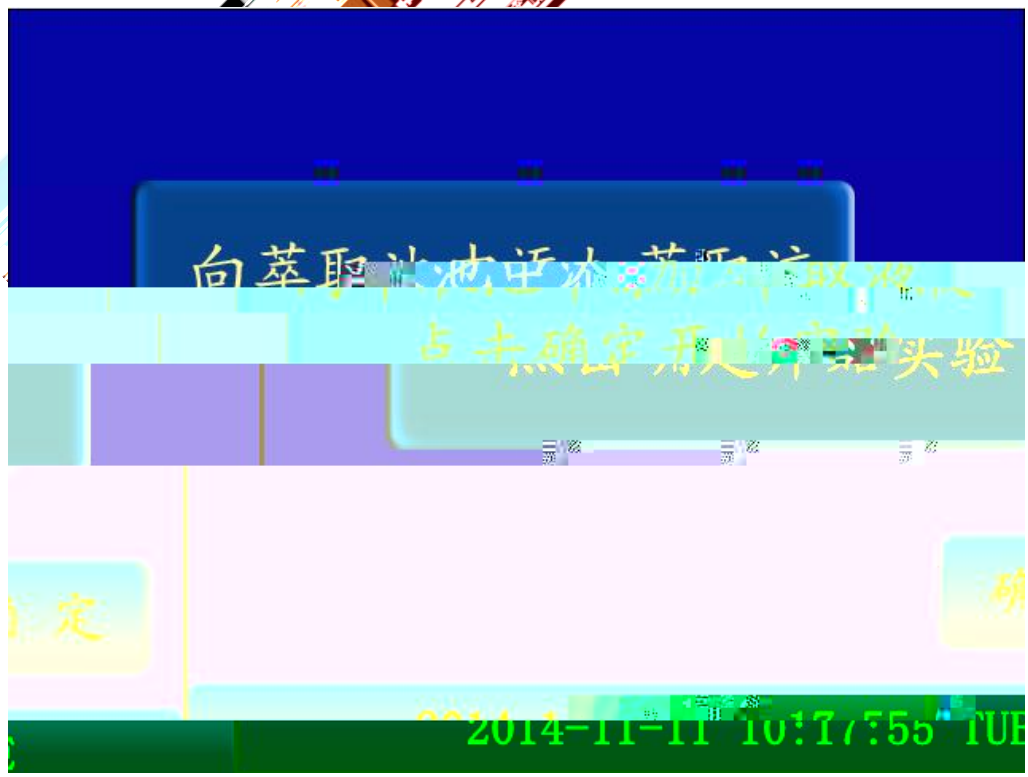


“ ”



正在中和萃取液
请稍候...

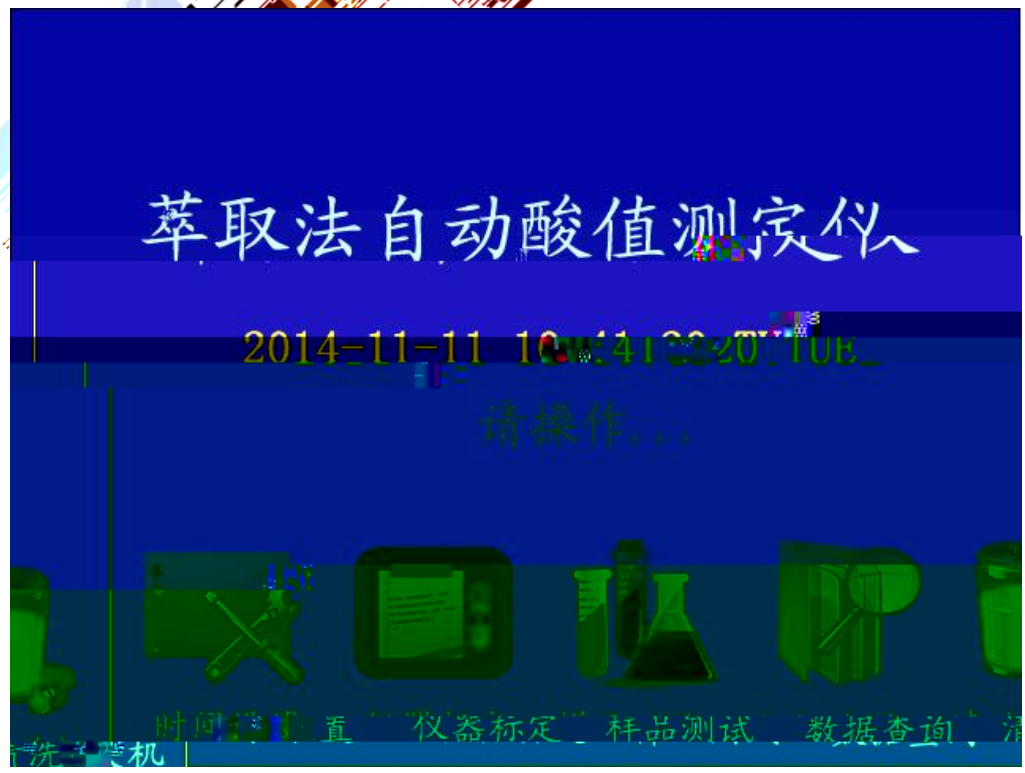
萃取法自动酸值测定仪
2014-11-11
请操作...
仪器标定
样品测试
数据查询

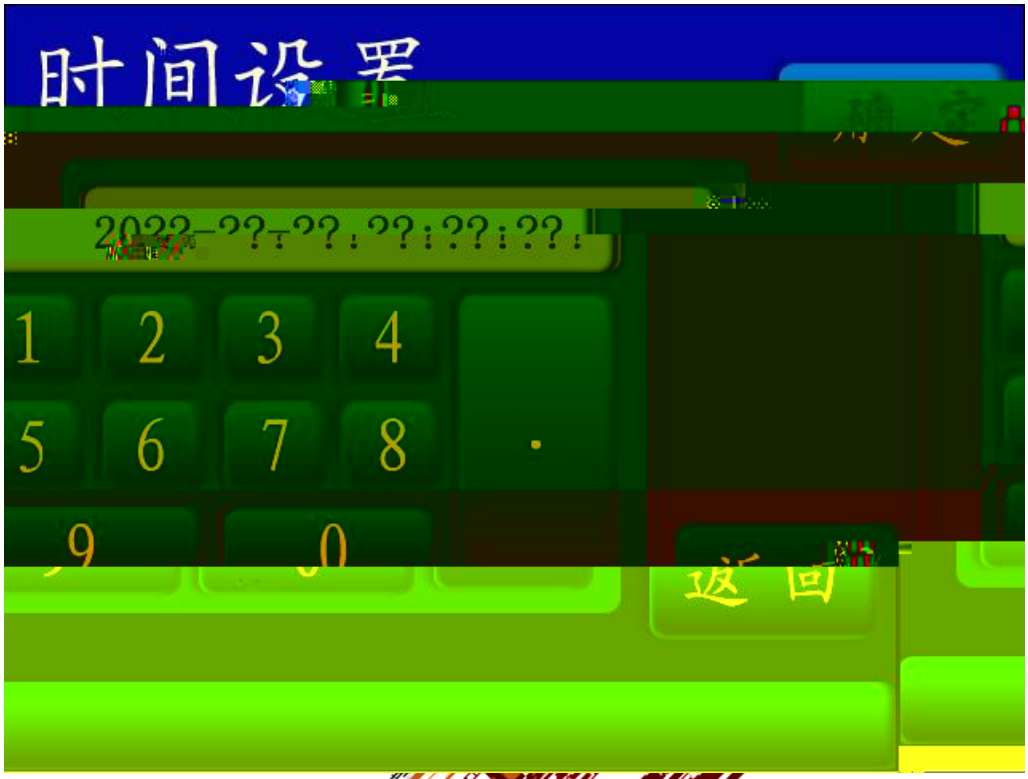
时间设置
2022-??-?? ??:??:??
1 2 3 4
5 6 7 8 .
9 0
返回

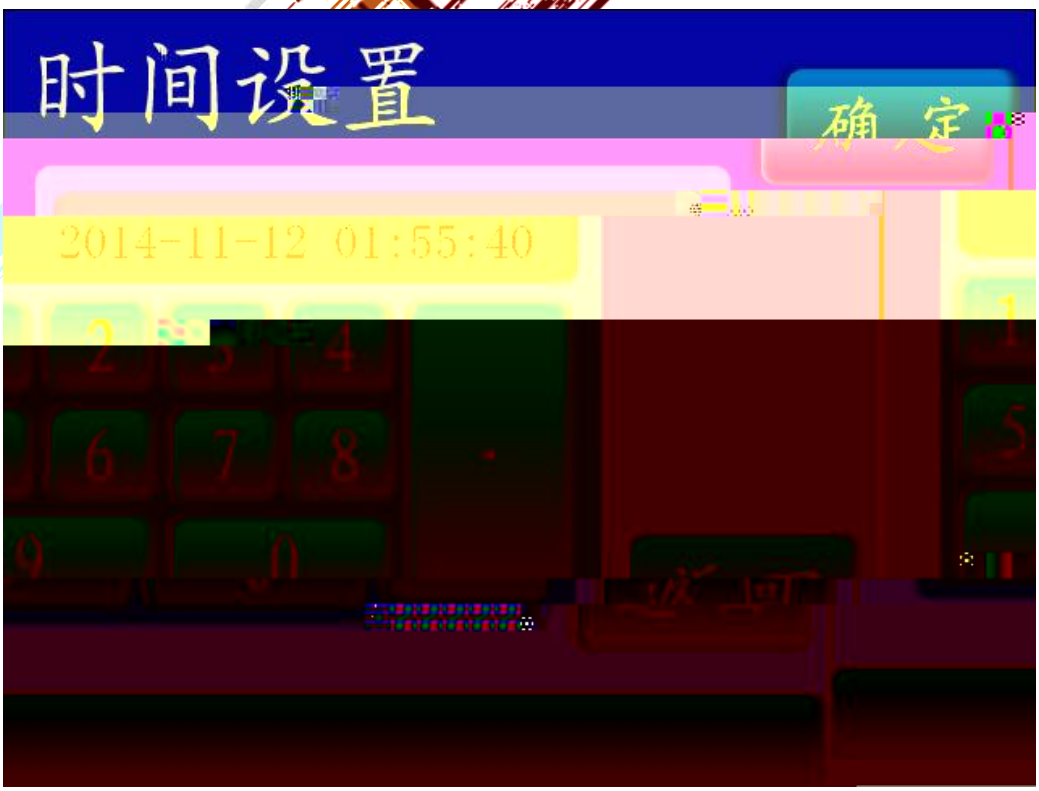

时间设置
确定
2014-11-12 01:55:40

“ ”

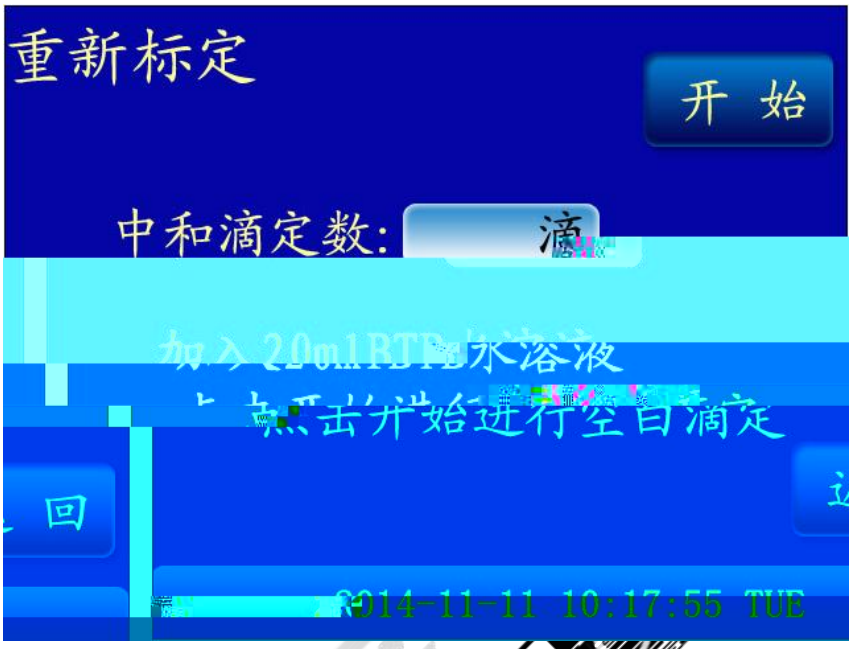


BTB 20ml BTB

BTB

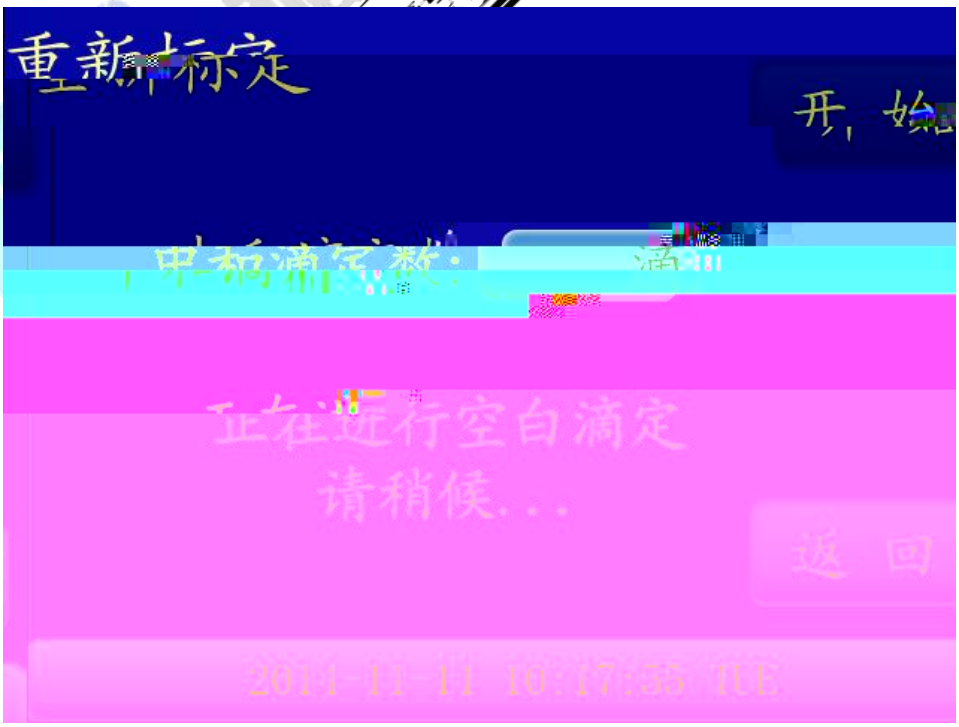

100 0.15mol/L


)

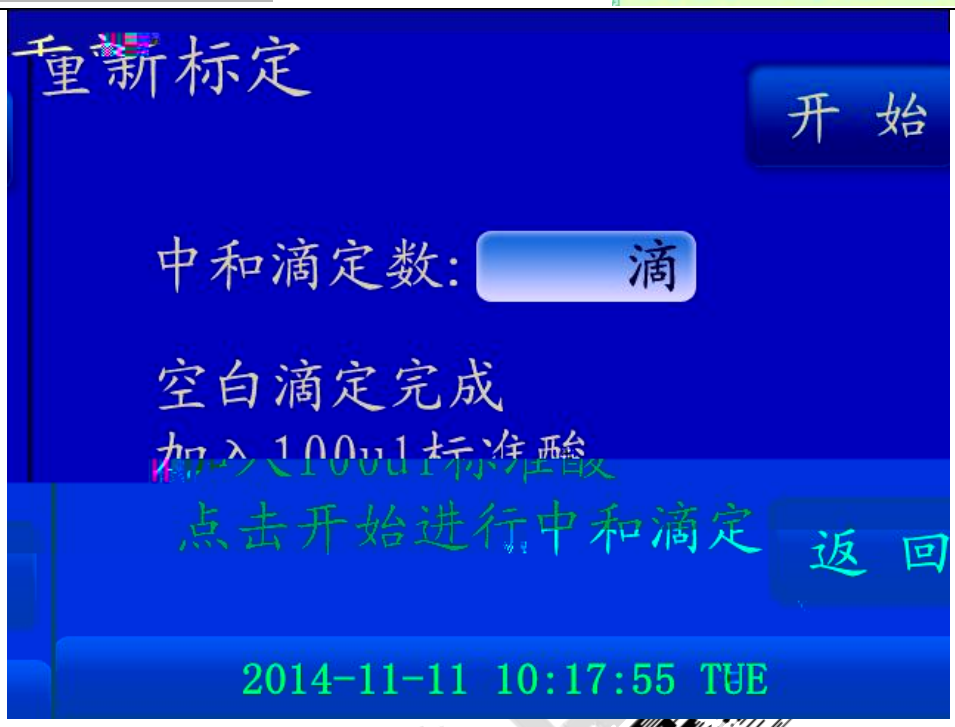
重新标定
开 始
中和滴定数: 滴
空白滴定完成
加入100ul标准酸
点击开始进行中和滴定
返 回
2014-11-11 10:17:55 TUE



重新标定
开 始
正在进行中和滴定
2014-11-11 10:17:55 TUE



重新标定
返 回



仪器标定
重新标定
10 滴
请选择
2014-11-11 10:17:55 TUE

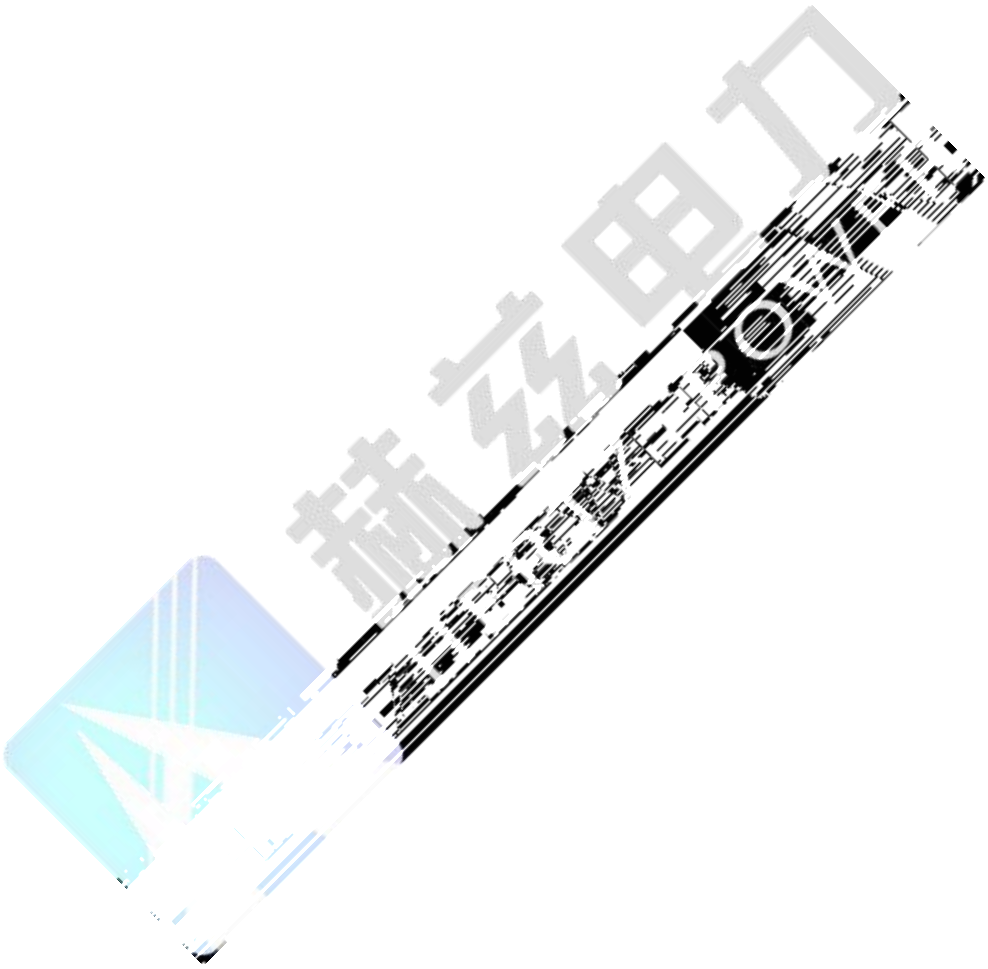


0.579
1 2 3 4
5 6 7 8
9 0
.
删 除
返 回

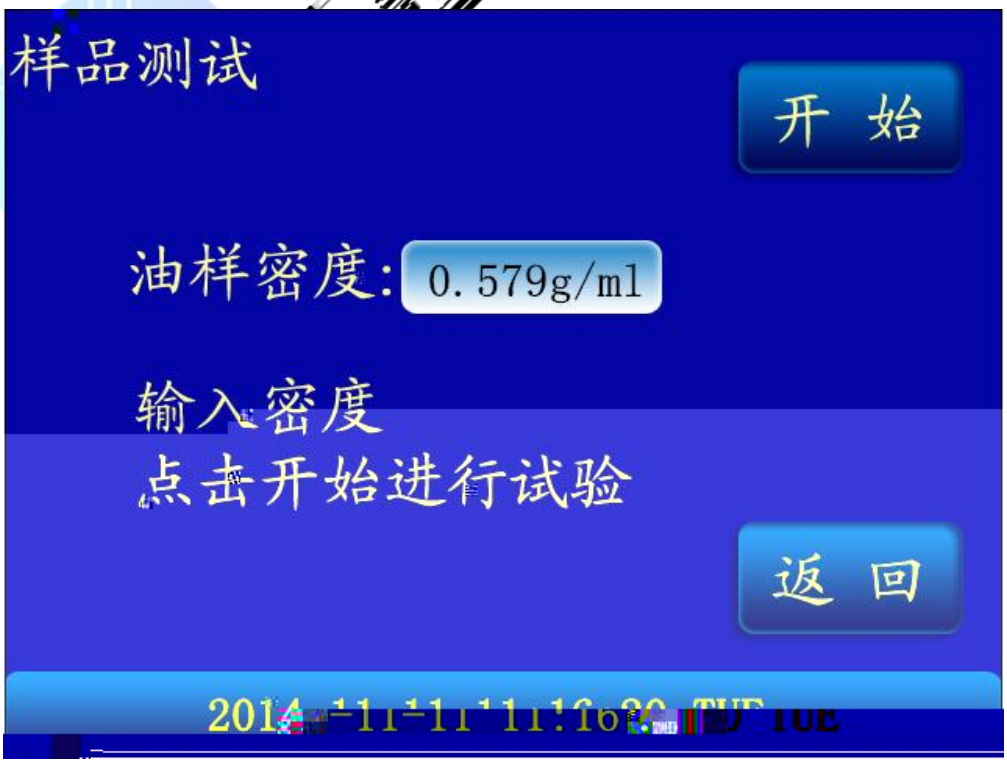

样品测试
开 始
油样密度:
0.579g/ml
输入密度
点击开始进行试验
返 回

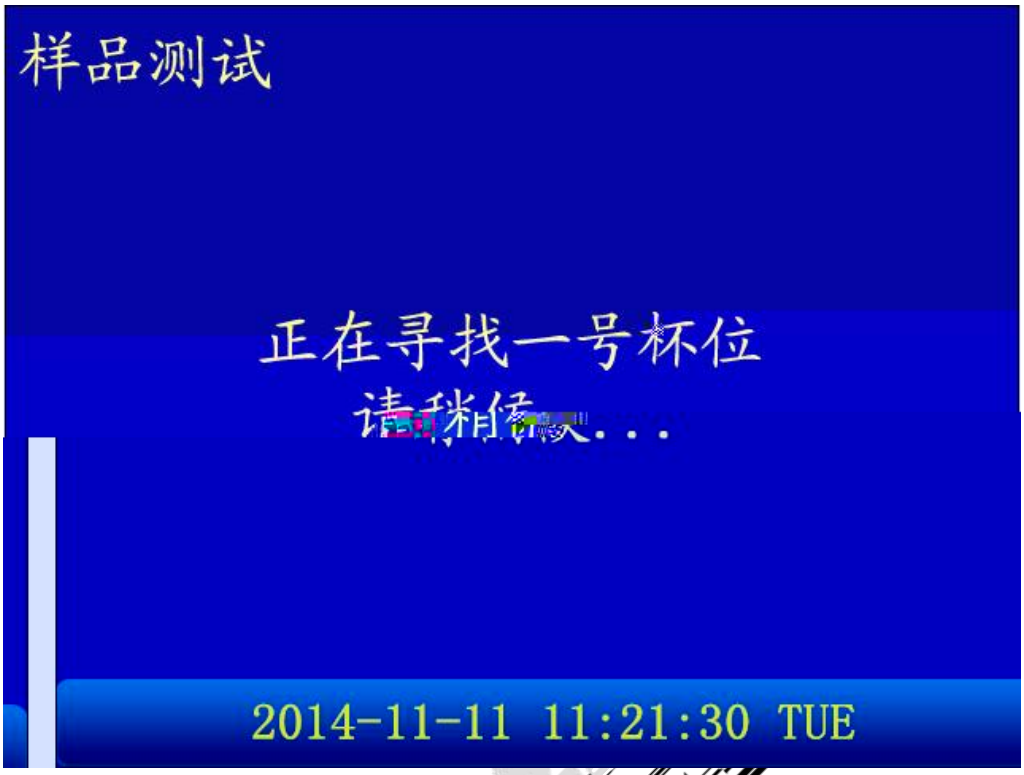
样品测试
正在寻找一号杯位
2014-11-11 11:21:30 TUE

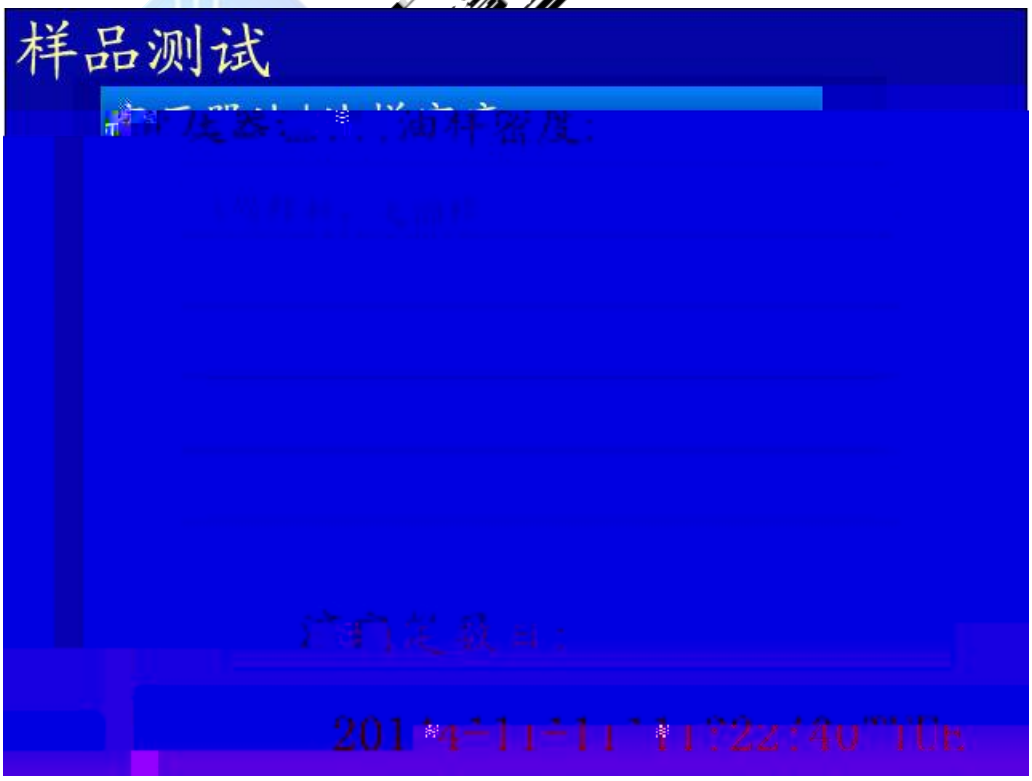

样品测试


10ml



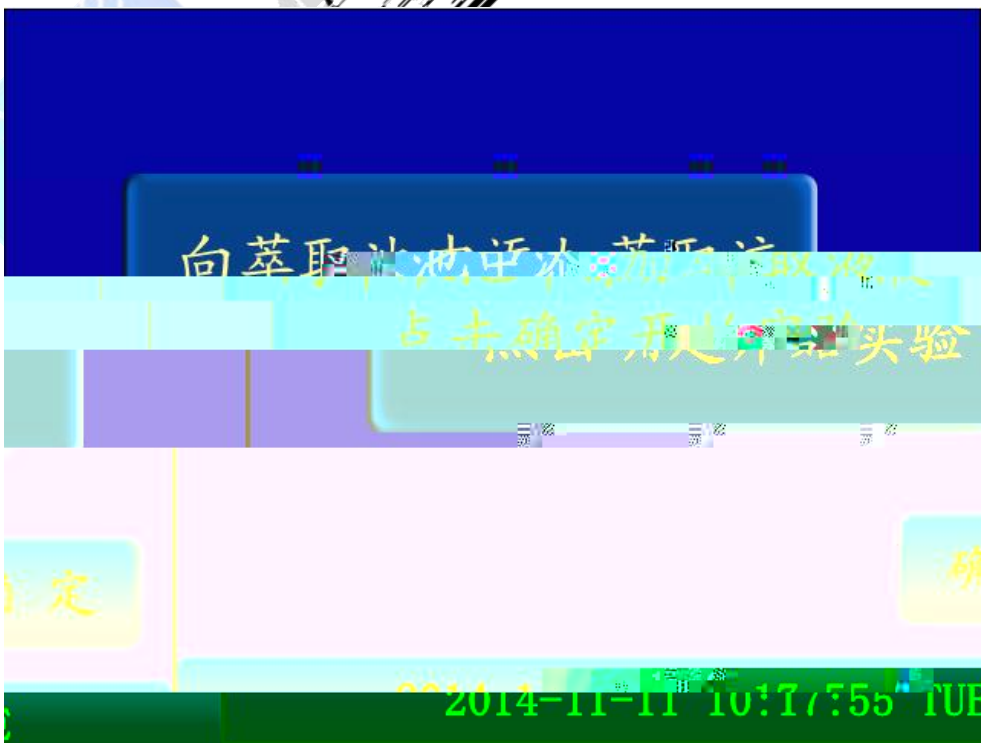

10ml

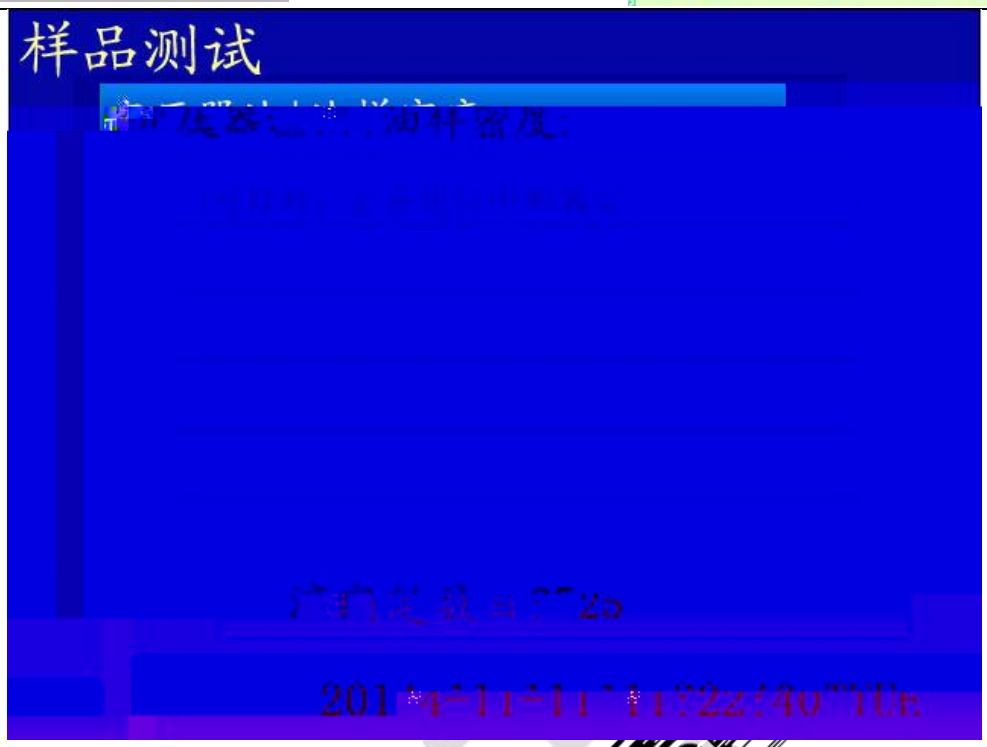
样品测试

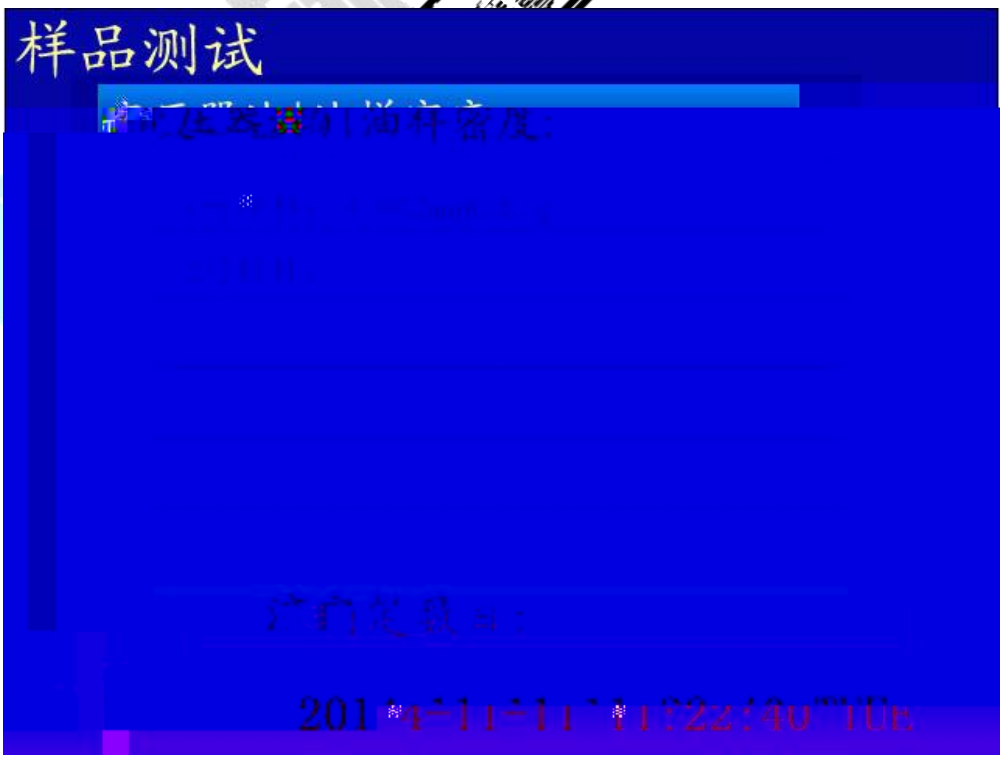
样品测试



样品测试
2号样杯： 无油样
3号油杯： 无油样
滴定数目：
2014-11-11 11:22:40 TUE


数据查询
无数据


数据查询
1 | 20
4号油杯： 0.002g/ml
5号油杯： 0.002g/ml
6号油杯： 无油样
删 除
返 回

“ ” “ ” “ ”


“ ” “ ”

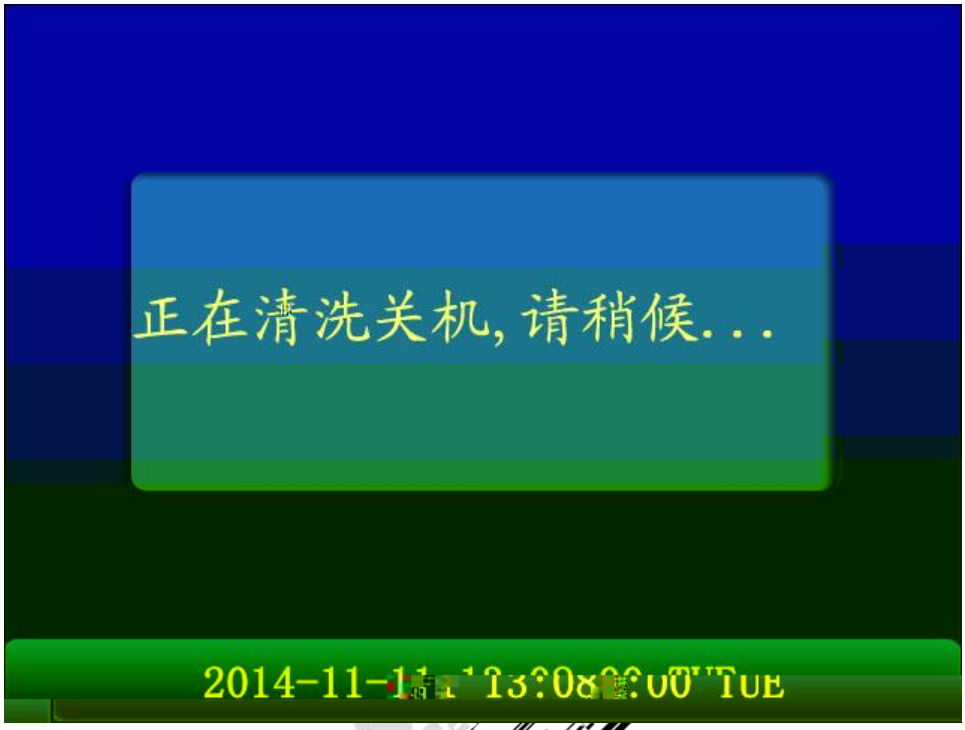
正在清洗关机，请稍候...

1

2

3 —KOH

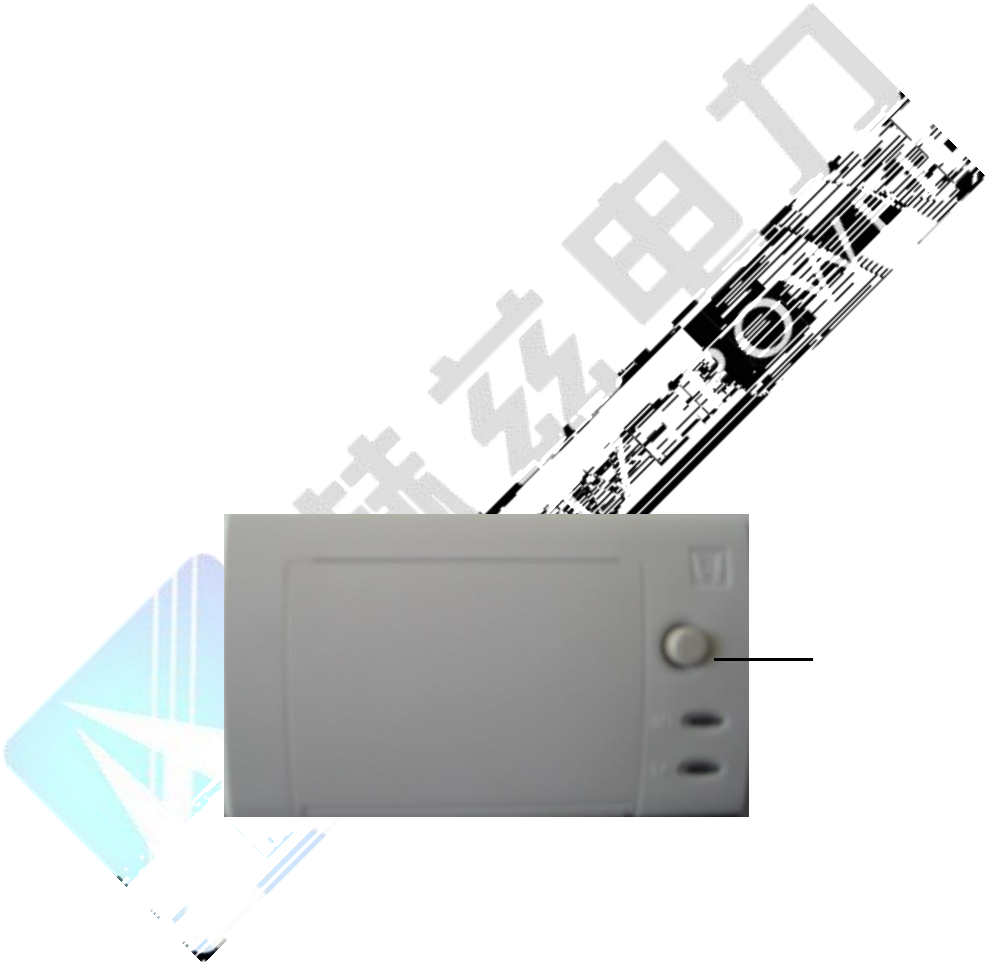

“ LF ” 5 “ LF ”

SEL

1

2

3

4

5 ,´

“

” “ ”

Ô 9 2 ‘ <¶$G° Y9 Đh ,DypÑÔ Ï ñ0 S

1 ˝ ≈

2 Ã Fs(Ÿ ñ0 ´Wi"f Ã

2

3 2 24

027-83267669

24 48

4

